



14 panel CLIA-waived drug test cup.

Branan ToxCup



Urine Drug Screening Device.

**CLIA
Waived**

For in vitro diagnostic use only.





How it works.

PRODUCT TRAINING AND CERTIFICATION AVAILABLE ONLINE.

Features & Benefits

- » CLIA-waived
- » Integrated test cup targets a wide range of illicit and prescription drugs
- » Simple procedure: collect specimen, tilt cup and read results in five minutes
- » Buprenorphine screening combined with low-level opiate cutoff makes it ideal for clinical settings
- » Self-contained cup is ideal for sending presumptive positive specimens for confirmation
- » Long shelf-life

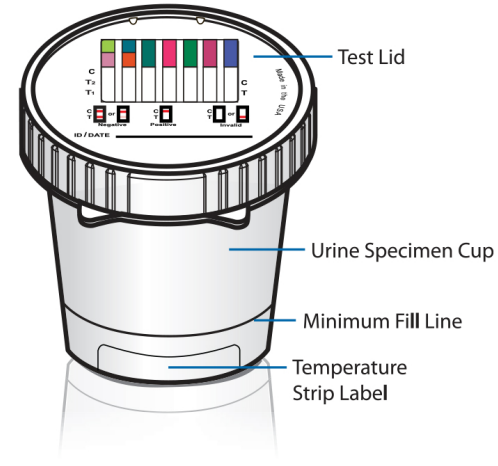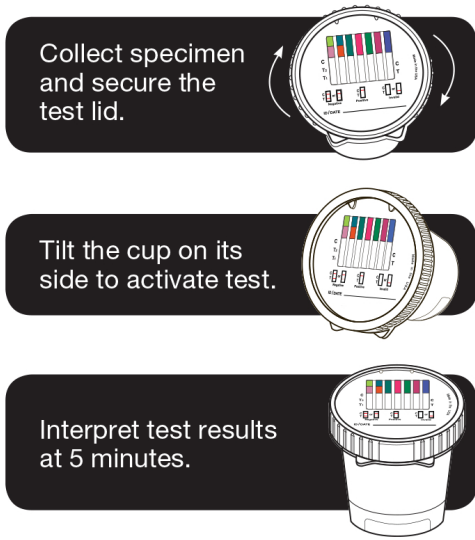

IMPORTANT PROCEDURE NOTE: Refer to product insert for complete instructions, limitations, and warnings. For in vitro diagnostic use only. Urine drug screen assay's provide only a preliminary analytical test result. A more specific alternate chemical method must be used in order to obtain a confirmed analytical result. Gas chromatography/mass spectrometry (GC/MS) or liquid chromatography/tandem mass spectrometry (LC/MS/MS) are the preferred confirmatory methods. Clinical consideration and professional judgement should be applied to any drug of abuse test result, particularly when preliminary positive results are indicated.

Drug Test Configuration

Item #	iCup® Dx 14
DT-14	AMP500, BAR300, BUP10, BZO300, COC150, MAMP500, MDMA500, MTD300, OPI300, OXY100, PCP25, PPX300, THC50, TCA1000

Distributed By:

U.S. SCREENING SOURCE

P.O. Box 676
 Pewee Valley, Kentucky 40056
 1 (866) 323-7336 – Toll Free
 (502) 412-4239 – Phone
 Fax Orders To: (502) 327-9037
 Email: customerservice@usscreeningsource.com
www.usscreeningsource.com